# Aprilla

## **AFM 2802**

Blackhead Suction Device

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
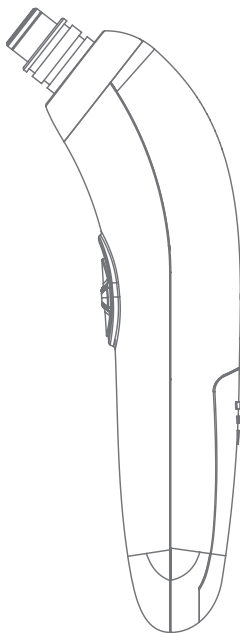
**EN FR DE IT ES PT**

## SAFETY INSTRUCTIONS

- By ignoring the safety instructions the manufacturer cannot be held responsible for the damage.
- This appliance shall not be used by children from 0 to 8 years old. This appliance can be used by children aged from 8 years old and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Keep the appliance out of reach of children aged less than 8 years old. Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised.
- **Batteries**
  - Do not expose the battery to high temperatures or direct sunlight. Never throw batteries into the fire. There is a danger of explosion!
  - Keep batteries away from children. They are not a toy!
  - Do not open the batteries by force.
  - Avoid contact with metallic objects.(Rings, nails, screws et cetera) there is a danger of shortcircuiting!
  - As a result of a shortcircuit batteries may heat up considerably or even catch fire. this may result in burns.
  - For your safety the battery poles should be covered with adhesive strips during transport.
  - If a battery leaks, do not rub your eyes with your hands. If in contact with the liquid, the hands should be washed thoroughly and eyes rinsed out with clear water, if symptoms persist please consult a doctor.

## PARTS DESCRIPTION

1. Small round hole head
2. Big round hole head
3. Oval head
4. On/Off switch
5. Battery cover
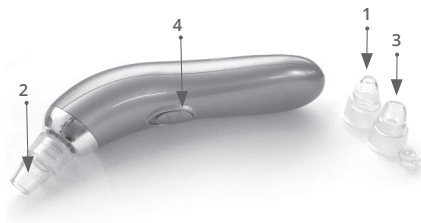
## BEFORE THE FIRST USE

### Battery installation

- Remove the battery cover by moving forward the little tab at the bottom of the cover and lifting it.
  - Insert 2 AA batteries in the battery compartment
- NOTE: MAKE SURE THE BATTERIES ARE PLACED IN THE RIGHT WAY LIKE INDICATED ON THE BATTERY COMPARTMENT.**
- Put the cover of the batteries compartment back on the device fitting first the cover top tab and then the bottom one. Press the cover until it is firmly closed.

## USE

- For best results the pores of the skin should be enlarged. For this reason it is a good idea to perform your skin care routine when taking your bath as heat from hot water will enlarge your pores.

- For best results, use the appliance two or three times a week as part of your skin care routine.

1. Remove any make-up from your face, leaving just bare skin. Never use the appliance on skin coated with cleansing lotion or facial soap.

2. Attach the Big round hole head to the head of the device, making sure that it is facing outward.

3. Move the appliance slowly around the tip of the nose.

4. Now turn the Power Switch off and connect the small round hole head. Then turn it on again.

5. Move the Small round hole head over the nose in the way described by the arrows. (As explained and shown in the Fig.)

6. Rinse your face, and apply moisturizing lotion or any other skin conditioner you normally use.

- Keep your face sufficiently moisturized around the tip and alas of the nose (Fig.)

1. Move the head all over the nose.

2. Limit the use of the Big round hole head over the whole nose to one minute or less.

3. Limit the use of the Small round hole head to two minutes or less. (Do not repeatedly move over the alas of the nostrils where the skin is thin.)

Use only the Small round hole head if you have concerns about skin irritation.

## CLEANING AND MAINTENANCE

The suction heads can be washed using water. Put some water in a bowl. (Do not use hot water and Do not use detergent) Wipe the body of the appliance and the heads with tissues or cotton buds, taking particular care to thoroughly wipe the insides of the suction head. Do not immerse the whole device in any liquid. Do not use detergent to avoid the risk of discolouration, deterioration, or cracking.

### GUIDELINES FOR PROTECTION OF THE ENVIRONMENT

This appliance should not be put into the domestic garbage at the end of its durability, but must be taken to a central point for the recycling of electric and electronic domestic appliances. This symbol on the appliance, instruction manual and packaging puts your attention to this important issue. The materials used in this appliance can be recycled. By recycling used domestic appliances you contribute to foment the protection of our environment. Ask your local authorities for information regarding the point of recollection.

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